


How to save and upload a segmentation

Here a simple guide on how to save and upload a segmentation

Save a segmentation


- 1) Apply the segmentation to the dataset
- 2) Use the command “Export segmentation” in the “Menu Options” under “File”
- 3) Select the “DICOM” option and save the segmentation
- 4) Move the segmentation file or files in a separate folder or directory. This folder has to contain only the segmentation files

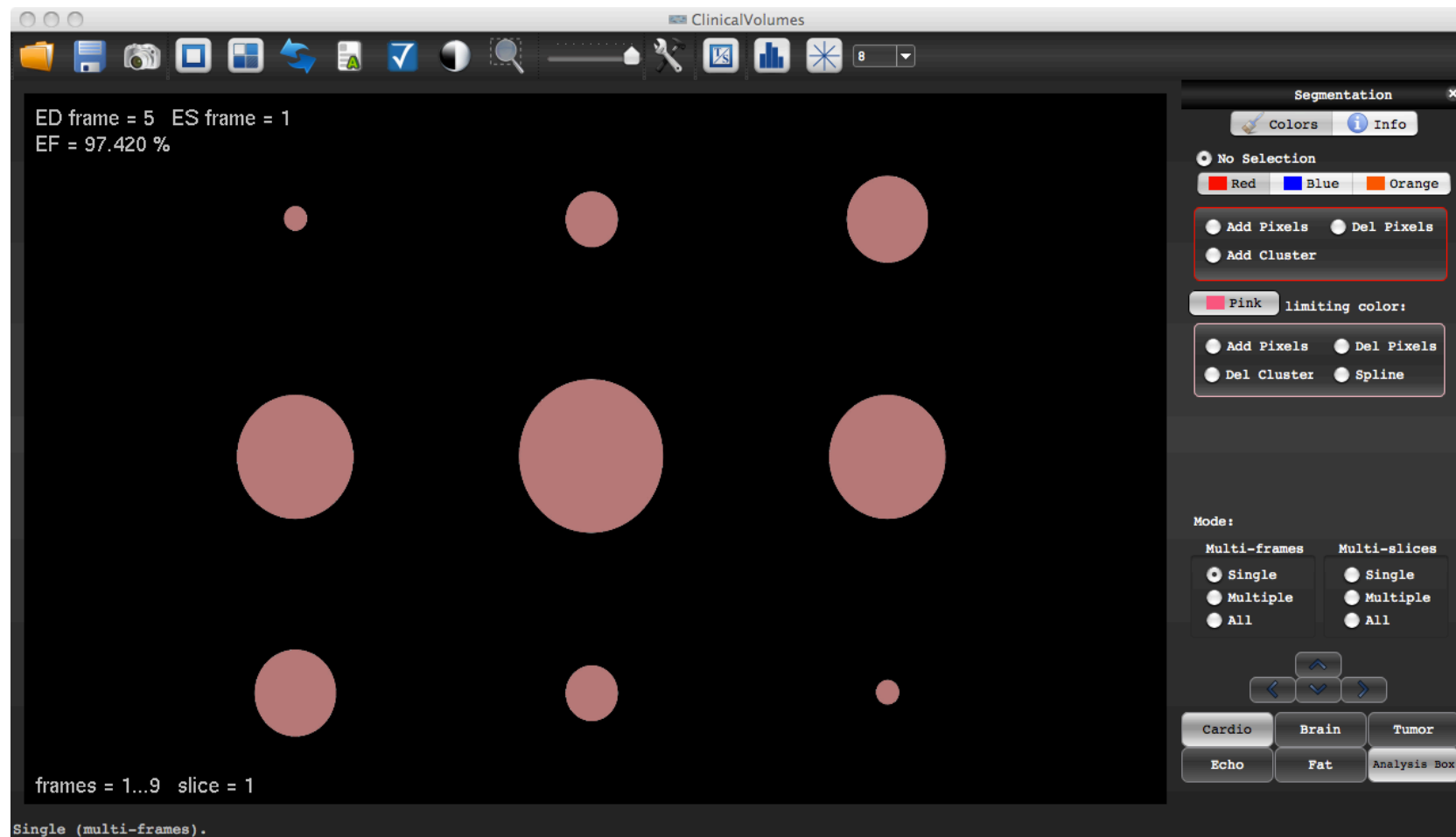

Load a segmentation

- 1) Open the original dataset (see “How to open a dataset”)
- 2) Use the command “Open segmentation” in the “Menu Options” under “File”
- 3) Select the folder or directory within which the segmentation file has been saved
- 4) Select “Choose”

Note: the segmentation files are dataset dependent. A segmentation file can not be uploaded on a different dataset rather than the one has been created from.

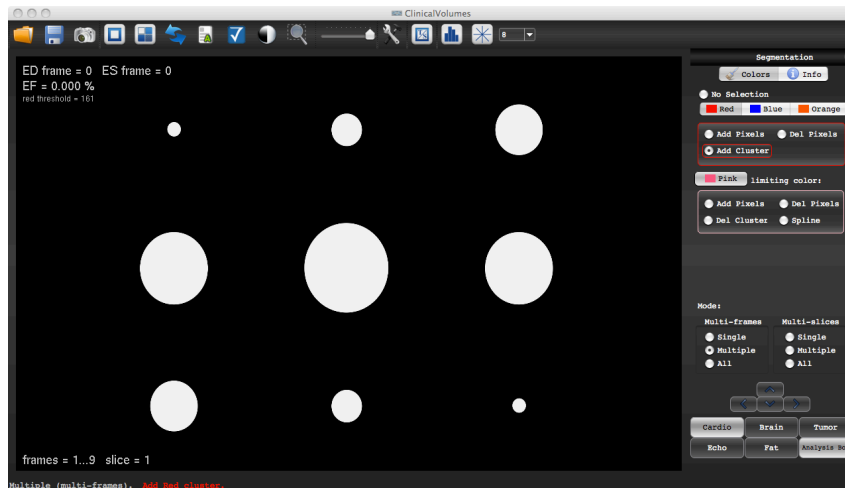

Open the original dataset (see “How to open a dataset”)

